

medi



lymphology
**measuring for
custom-made
inelastic circaid[®]
lower extremity**

medi. I feel better.

circaid® Custom Legging Form

☐ Exact Reorder Order

Number: _____



Date _____ Physicians Printed Signature _____

NPI _____ Physicians Signature _____

Customer Name _____

Customer No. _____ Purchase Order No. _____

Patient Name _____ Fax _____

Billing Address _____

Shipping Address _____

Measured by _____

Telephone _____

Order Date _____ Email _____

Diagnosis: ICD-10 - check all that apply

- ☐ I89.0 Lymphedema, not elsewhere classified
☐ I97.2 Postmastectomy lymphedema syndrome
☐ I97.89 Other postprocedural complications and disorders of the circulatory system, not elsewhere classified
☐ Q82.0 Hereditary lymphedema

Shipping method: ☐ Standard (max. 5 days after complete order is received) ☐ Second Day (extra charge) ☐ Next Day (extra charge)

Credit Card Info _____

juxtafit® premium					
	Qty Left	Qty Right	Lateral rise (oblique)	Add pull tabs	Add Ez-on panel
Lower leg			☐ No rise (default) ☐ 5cm ☐ 10cm		
Knee only					
Upper leg			☐ Yes (default) ☐ No		
Upper leg w/knee			☐ Yes (default) ☐ No		
Whole leg (one piece)			☐ Yes (default) ☐ No		
Custom foot					

Additional Options		
☐ pac band™ compression anklets included ☐ single band afw™ ☐ juxtafit premium interlocking afw ☐ juxtafit premium afw ☐ customizable interlocking afw	Extra Pair Undersleeves (open ended): Lower Leg _____ Whole Leg _____ Beige _____ Silver _____ Extra Pair Undersocks (close ended): Lower Leg _____ Whole Leg _____ Beige _____ Silver _____ Cotton Terry _____	
Garment color ☐ beige (default) ☐ black	Compressive Undersocks 15-25 mmHg 25-35 mmHg Foot length Small _____ Small _____ ☐ Standard Large _____ Large _____ ☐ Long Foot Cover up ☐ black ☐ beige	

graduate™								
	Qty Left	Qty Right	Boot style	Foam lateral rise (oblique)	Band locks	Hard sole boot	Foam liner closure	Foam liner color options
Lower leg with boot			☐ attached (default) ☐ separate					Interior ☐ beige *
Lower leg without boot								Exterior ☐ black* ☐ beige
Whole leg with boot			☐ attached (default) ☐ separate					* default
Whole leg without boot								
Boot only								
			Foam pad accessory: _____ cm x _____ cm (max. 20cm x 20cm)					

Notes on taking measurements for custom-made circaid inelastic products

• Measurements for compression garments should not be taken until the best possible decongestion effort has been achieved. Circumference and length measurements are taken in a distal to proximal sequence. Measurements should be taken with the leg extended straight, such as with a patient lying flat on a table or standing.

• It is essential to mark the measuring points on the leg so that the circumference and length measurements are taken at the same point.


• The circumference measurements are taken without any tension as the products are adjustable and will accommodate some changes in size. Measurements are to be taken with skin measurements. Skin measurements should be taken loosely without tension.

Patient Name: _____

Date: _____



Foot measurement Circumferences

Left in cm.		Right in cm.
	c _{mall} .	
	c _y	
	c _{A1}	
	c _A	

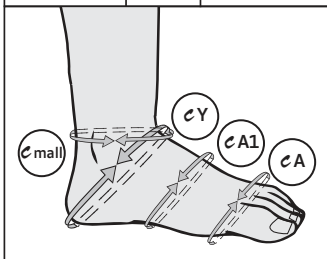

Foot measurement Lengths






Left in cm.		Right in cm.
	ℓ _{P2}	
	ℓ _{floor-mall.}	
	ℓ _{A-Y}	
	ℓ _A	

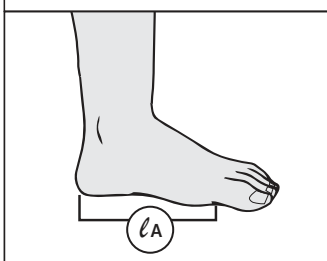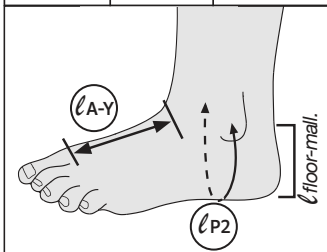

Please take measurements without tension!

Leg measurement Circumferences

Left in cm.		Right in cm.
	85	
	80	
	75	
	70	
	65	
	60	
	55	
	50	
	45	
	40	
	35	
	30	
	25	
	20	
	15	
	10	
	05	
	c _{mall}	
	c _{E*}	

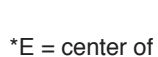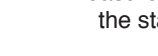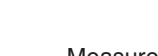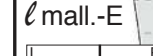ℓ_{mall.-K}


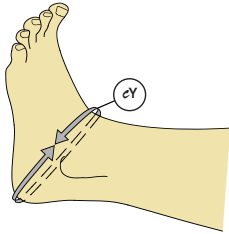
Measurements must be every 5cm from the starting point at the malleolus.

*E = center of patella


Notes:

Foot & Ankle Measurements


All measurements are taken in a straight line.




1. Circumference (c) Y - Measurements taken over the heel and arch in dorsal flexion with tendons taut; take physical dimensions without tension; Take soft tissue projections into account.

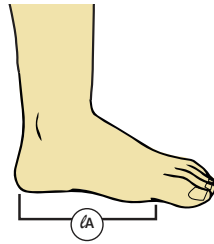


2. Length (l) A - Length from the base of toes to the back of the heel. (medial length)

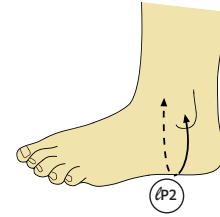

3. Length (l) P2 - Length from malleolus to malleolus; pass under heel.

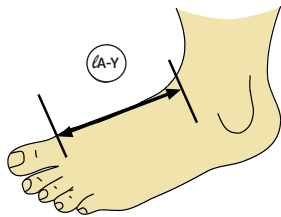


4. Length (l) A-Y - Length from ball of foot to where foot meets leg.

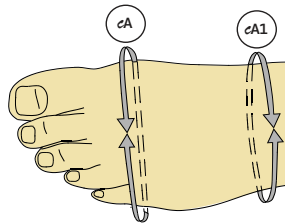

5. Circumference (c) A - Circumference at ball of foot.

Circumference (c) A1 - Circumference at middle of foot.

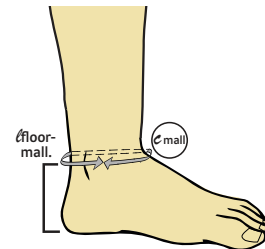


6. Circumference (c) mall. - Circumference at malleolus.


Length (l) floor-mall. - Length from floor to malleolus.

Leg Measurements


All measurements are taken in a straight line. Measurements start from malleolus, not floor. If shelves, folds or trouble areas are present please draw shape of patient's leg darker on sample or separate sheet. Pictures and contoured measurements are also helpful for challenging limb shapes.

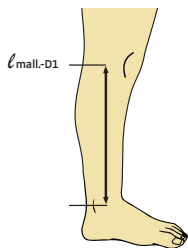


7. Length (l) mall-D1 - Medial length from malleolus to knee crease (where product will typically begin & end)

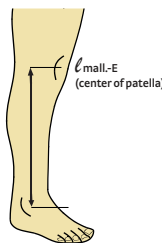

8. Length (l) mall-E - Lateral length from malleolus to center of patella (required for all except lower leg).

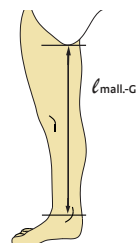

9. Length (l) mall-G - Medial length from malleolus to groin, far enough below pubic area that garment will not cause discomfort (required for all except lower leg). (where product will typically begin & end)

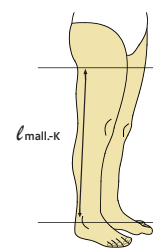

10. Length (l) mall-K - Lateral length from malleolus to gluteal fold (required for all except lower leg).